
LabActor

LabActor main control pannel

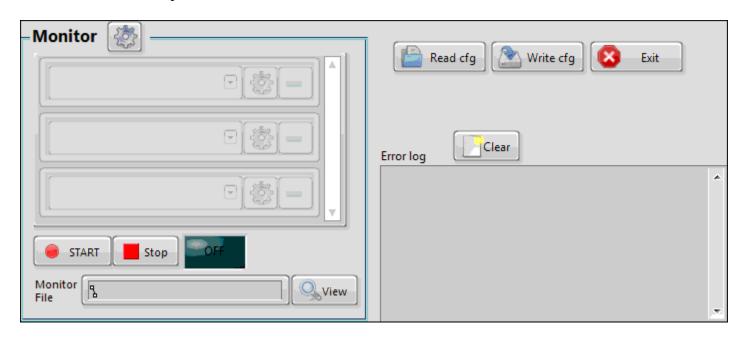

[abc] Monitor

List of active measurement modules

TFI Start monitor

Start measurement

III Stop monitor

Stop measurement

TEI View XY Plot

View Monitor data file

TFI Read cfg

Read application configuration from file

TFI Write cfg

Write application configuration to file

III Monitor configuration

Open Monitor configuration

TFI Exit

Exit application

TFI Clear log

Clear error log

Monitor Running

Measurement state


Monitor File

Latest Monitor data file

Error log

Monitor

Monitor configuration

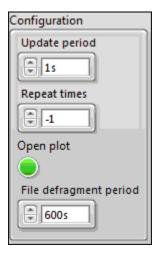

Configuration

Monitor configuration

DBLI Update period

Request data collection this often

1321 Repeat times

Complete this many measurements then stop (or keep going if <0)

TFI Open plot


Open new Monitor data plot on start

File defragment period

Defragment Monitor data file this often when running

Spectrum Recorder

Spectrum Recorder module sweeps frequency and records complex response using Red Pitaya

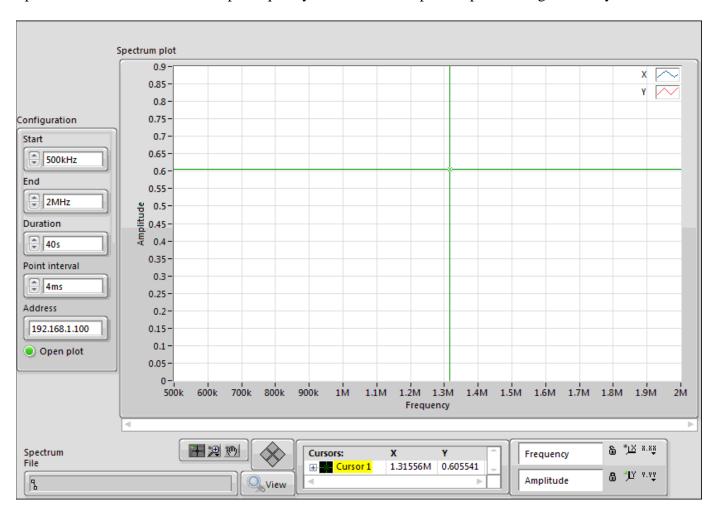

Essis Configuration

Spectrum recorder configuration

DBLI Start

Spectrum starting frequency

DBLI End

Spectrum ending frequency

Duration

Spectrum duration

DBLI Point interval

Interval between frequency points

Address

Network address of Red Pitaya

TFI Open plot

Open new spectrum plot window at the start of measurements

TFI Regions of interest

[DBL] Region centers

DBL Numeric

DBLI Regions width

TIP View XY Plot

View latest spectrum file

Spectrum plot


Latest spectrum

Spectrum File

Latest spectrum file

Line Detector

Line Detector module finds amplitude peaks in spectrum file and fits complex Lorentzian around peaks

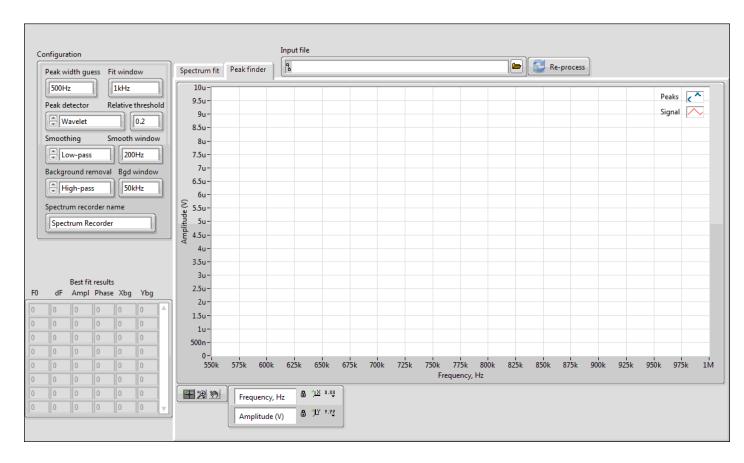

Input file

Spectrum file to process (auto-updates when monitor is running)

Configuration

Line Detector configuration

DBL Peak width guess

Guess width for peak finder

Background removal

Background removal type for peak finder

Smoothing

Smoothing type for peak finder

DBL Smooth window

Smoothing amount for peak finder

DBLI Relative threshold

Smallest peak magnitude relative to spectrum maximum

Peak detector

Peak finder type

Bgd window

Background removal width for peak finder

Spectrum recorder name

Name of the spectrum recorder module to watch

IBLI Fit window

Width of complex Lorentzian fit window around the peak

TFI Re-process

Re-process spectrum selected in the Input file window

Spectrum plot

Spectrum data with complex Lorentzian fit

Peak Plot

Filtered and smoothed magnitude with peaks found

[DBL] Best fit results

Latest complex Lorentzian fit results

Numeric